
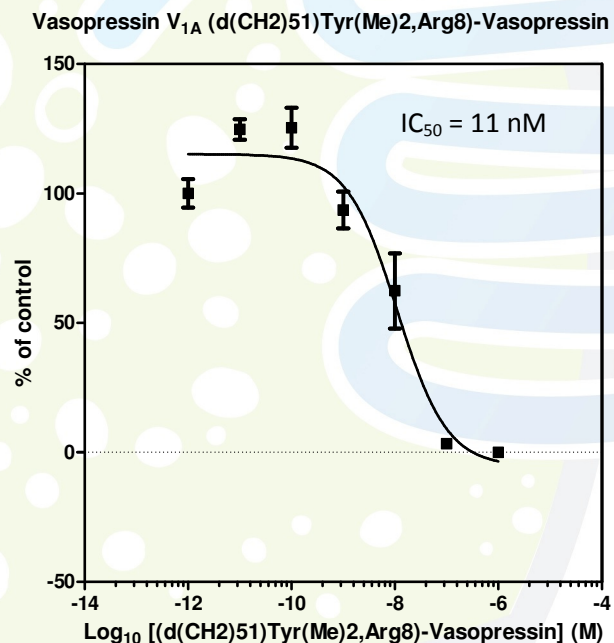


## GPCR Profiling

Target:	Vasopressin V1a
Available Assays:	Radioligand Binding & Cell-Based
Cell-Based Agonist:	Arg8-Vasopressin
Cell-Based Antagonist:	(d(CH2)51Tyr(Me)2,Arg8)-Vasopressin

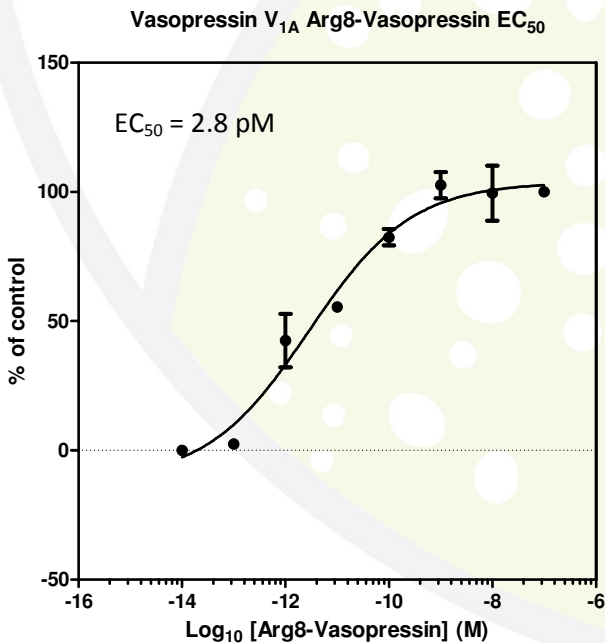

GPCR-expressing cells were plated in black, clear-bottomed 96 well plates at a density of 25,000 cells per well in 100µl volume and incubated overnight. Next day, medium was removed and 25 µl HBSS buffer (HEPES buffered Hank's balanced salt solution) was added to the wells. 10 µl calcium 5 dye (prepared according to the manufacturer's instructions in HBSS) was added and plates incubated at 37°C for 1 hour. If appropriate, antagonists were added and incubated for 5 min at RT. The plates were then placed in a fluorescence reader whereupon agonist was added and fluorescence monitored.